

CUSTOMER: _____

REFERENCE: _____

GETINGE

GETINGE GROUP

4003 FLOOR LOADING CARTS

PRODUCT SPECIFICATION

PRODUCT

The Getinge Floor Loading Sterilizer Cart is an all-purpose cart designed for processing, transportation and storage. The cart is constructed to resist the corrosive action of high temperature steam sterilization. The cart is specifically designed for use with Getinge Floor Loading Sterilizers and Washers. Since the cart is mobile, it saves space that would otherwise be dedicated for conventional transfer carriages.

SIZE

Length	Width	Height
36 ³ / ₄ "	24 ³ / ₄ "	57 ⁵ / ₈ "
934mm	629mm	1463mm

MODELS

- ☐ 4003: Includes 2 shelves/rods and swivel casters
- ☐ 4003RC: Includes 2 shelves/rods and fixed casters
- ☐ WCY61301604959: Includes no adjustable shelves or retaining rods, does include four swivel casters

QUALITY STATEMENT

Confidence in the Getinge Group is the most important quality criterion. This must be the hallmark of all our external and internal commitments, activities and products. Products and services supplied by Getinge must conform to the agreed terms and expectations to ensure recommendations for further business. The achievement of these quality goals is the basis for a continued competitive and successful enterprise.

CONSTRUCTION

The entire cart is constructed of 300 Series stainless steel. Frame, tubes and rungs are constructed to enable the cart to carry loads of up to 500 pounds (226 Kg). The Floor Loading Sterilizer has interior chamber rollers that are designed to engage the lower frame of the cart for easy sliding into and out of the Floor Loading Sterilizer and to keep the wheels above the chamber floor. This applies to both pit mounted or floor mounted sterilizers. The wheels are constructed from Teflon® material, which is capable of withstanding repeated sterilization without flat spotting and require no lubrication. There are four swivel casters and two kick-down locks for tracking. A locking bar is provided at one end to enable the unloading of two carts. When the Floor Loading Sterilizer is floor mounted the Model 4003RC Cart is designed to interface with the 4004 Transfer Carriage. The Model 4003 Cart will interface with the hydraulic lift fixed permanently in front of the floor mounted sterilizer door. Bumpers are provided at each corner. Four retaining rods are provided.




SHELVING

The cart is supplied with a fixed base shelf and two adjustable shelves. Shelves are constructed with heavy frames and wire rods. End rods are spaced at even intervals for easy shelf adjustments.

ACCESSORIES

- ☐ Extra Shelves: Cat #CSY61301600641 Quantity _____
- ☐ Retaining Rods: Cat #FSY61301604960 Quantity _____

FLOOR LOADING CARTS

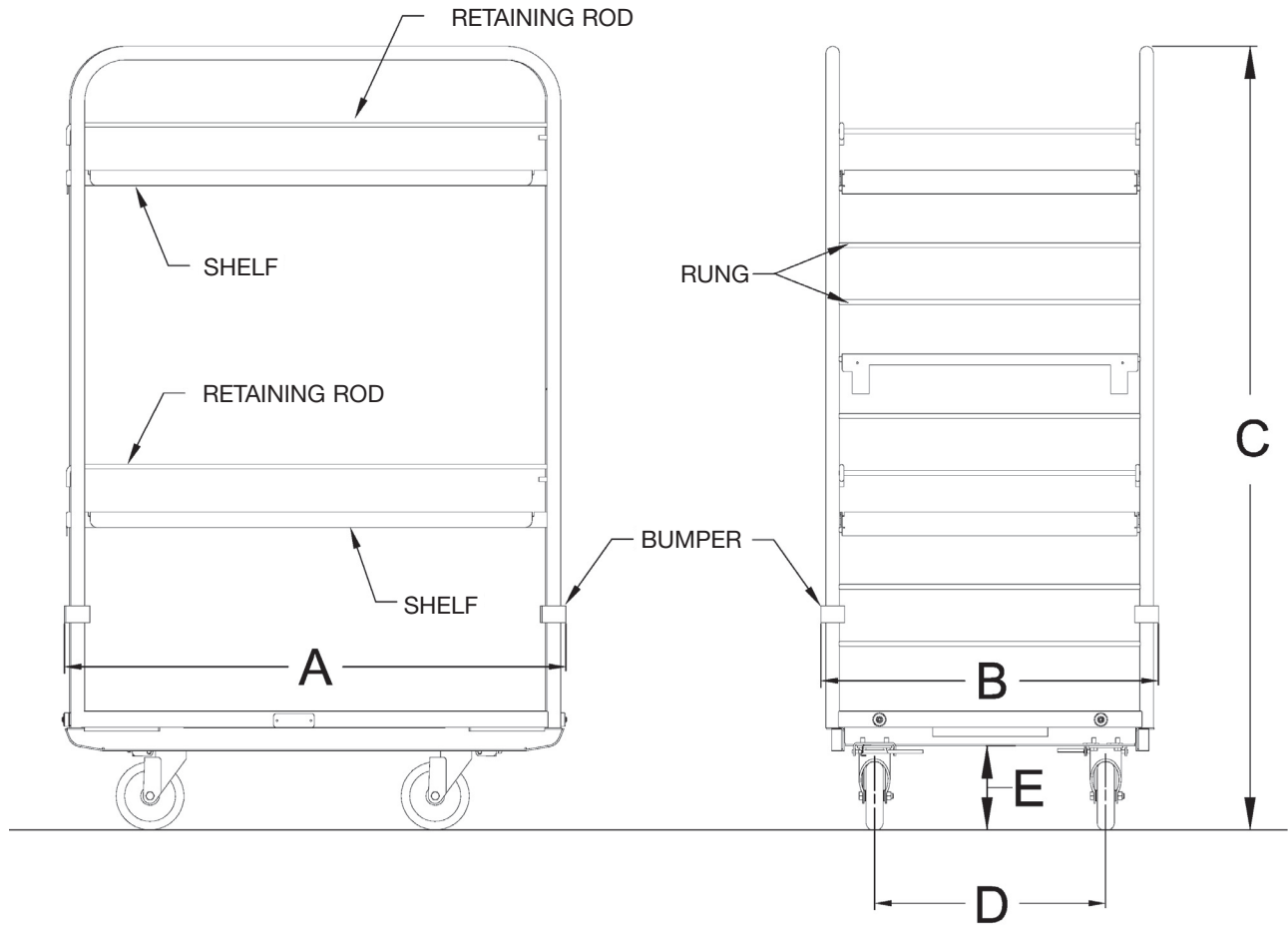

Model	A	B	C	D	E	Rung Spacing
4003	36 ³ / ₄ " (934 mm)	24 ³ / ₄ " (629 mm)	57 ⁵ / ₈ " (1463 mm)	17" (432 mm)	6 ¹ / ₄ " (159 mm)	4 ³ / ₁₆ " (107 mm)

Geringe USA, Inc
 1777 East Henrietta Road
 Rochester, New York 14623-3133
 USA
 Phone: (800) 475-9040
 Fax: (585) 272-5033
info@getingeusa.com
www.getinge.com

Geringe Canada Ltd.
 1575 South Gateway Road, Unit C
 Mississauga, Ontario L4W 5J1
 Canada
 Phone: (905) 629-8777
info@getingeusa.com
www.getinge.com

GETINGE
 GETINGE GROUP

www.getinge.com

GETINGE GROUP is a leading global provider of equipment and systems that contribute to quality enhancement and cost efficiency within healthcare and life sciences. We operate under the three brands of ArjoHuntleigh, GETINGE and MAQUET. **ArjoHuntleigh** focuses on patient handling and hygiene, disinfection, DVT prevention, medical beds, therapeutic surfaces and diagnostics. **GETINGE** provides solutions for infection control within Healthcare and contamination prevention within Life Sciences. **MAQUET** specializes in therapeutic applications, products, solutions and services for OR and ICU.